Medical Gas Maintenance Report

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Oxygen - Nitrous Oxide - Air - Vacuum - Liquid Oxygen Tank - Oxygen Generator

Cairo - Egypt : 21 Hafez Al - Abhar Street
02 / 4872003 02 / 26821616

01271332581 - 01032744799 - 01009084613 - 01015013648

Giza - Egypt : El Nozha Boulak El Dakrour

MEDICAL GAS MAINTENANCE CHECK LIST

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MANUFACTURER →

DATE OF MANUFACTURER →

	COMPRESSOR # 1	COMPRESSOR # 2	COMPRESSOR # 3
MODEL			
SERIAL NUMBER			
POWER NEW			
FLOW			
CONTROL PANEL			
OIL FILTER			
OIL LEVEL			
RUNNING HOURS			

DRYERS

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	DRYER # 1	DRYER # 2	DRYER # 3
MANUFACTURER			
SERIAL NUMBER			
DRAIN			
FLOW			
TYPE			
POWER			

TANK

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MANUFACTURER	
CAPACITY	
CLEAN	
DRAIN	
SAFETY VALVE	

FILTERS

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CYCLONIC FILTER	
MICRONIC FILTER	
SUBMICRONIC FILTER	
DUST FILTER	
AUTOMATIC DRAIN	

REGULATORS

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4 BAR	
7 BAR	

RECOMMENDATION

NOTES

IGAS ENG NAME:

CONT ENG NAME:

M.O.H ENG NAME:

IGAS ENG SIGNATURE:

CONT ENG SIGNATURE:

M.O.H ENG SIGNATURE:

MEDICAL GAS MAINTENANCE CHECK LIST

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MANUFACTURER

DATE OF MANUFACTURER

	PUMP # 1	PUMP # 2	PUMP # 3
MODEL			
SERIAL NUMBER			
POWER			
OIL FILTERS			
OIL LEVEL			
CONTROL PANEL			
FLOW			
RUNNING HOURS			
MAX PRESSURE			
MIN PRESSURE			

TANK	التانك
MANUFACTURER	
CAPACITY	
DRAIN	

FILTERS	الفلاتر
MANUFACTURER	
CLEAN	
TYPE	

RECOMMENDATION

NOTES

IGAS ENG NAME:
CONT ENG NAME:
M.O.H ENG NAME:

IGAS ENG SIGNATURE:
CONT ENG SIGNATURE:
M.O.H ENG SIGNATURE:

MEDICAL GAS MAINTENANCE CHECK LIST

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MANUFACTURER	SERIAL NO
TANK CAPACITY	DATE OF MANUFACTURING
LEAK TEST	GAUGE
OK ☐ DEFEL ☐	OK ☐ DEFEL ☐
VAPORIZER	REFILL TIME
OK ☐ DEFEL ☐	

RECOMMENDATION
IGAS ENG NAME:
CONT ENG NAME:
M.O.H ENG NAME:

NOTES
IGAS ENG SIGNATURE:
CONT ENG SIGNATURE:
M.O.H ENG SIGNATURE:

MEDICAL GAS MAINTENANCE CHECK LIST

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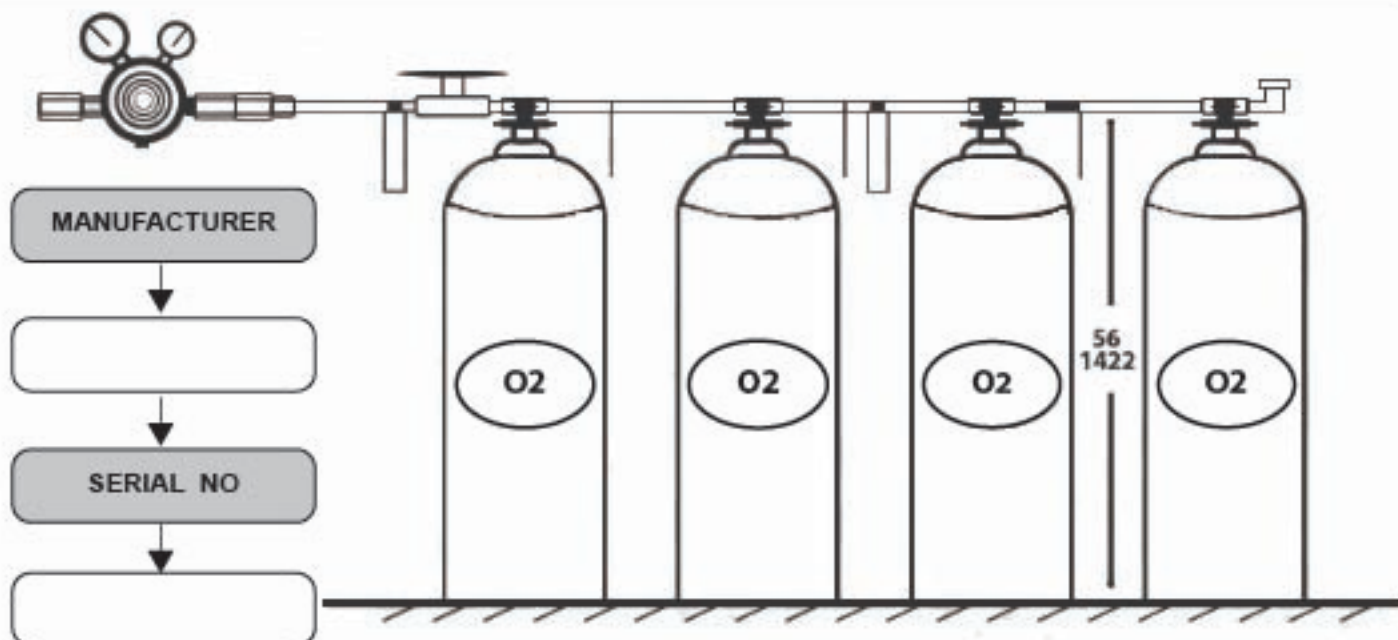
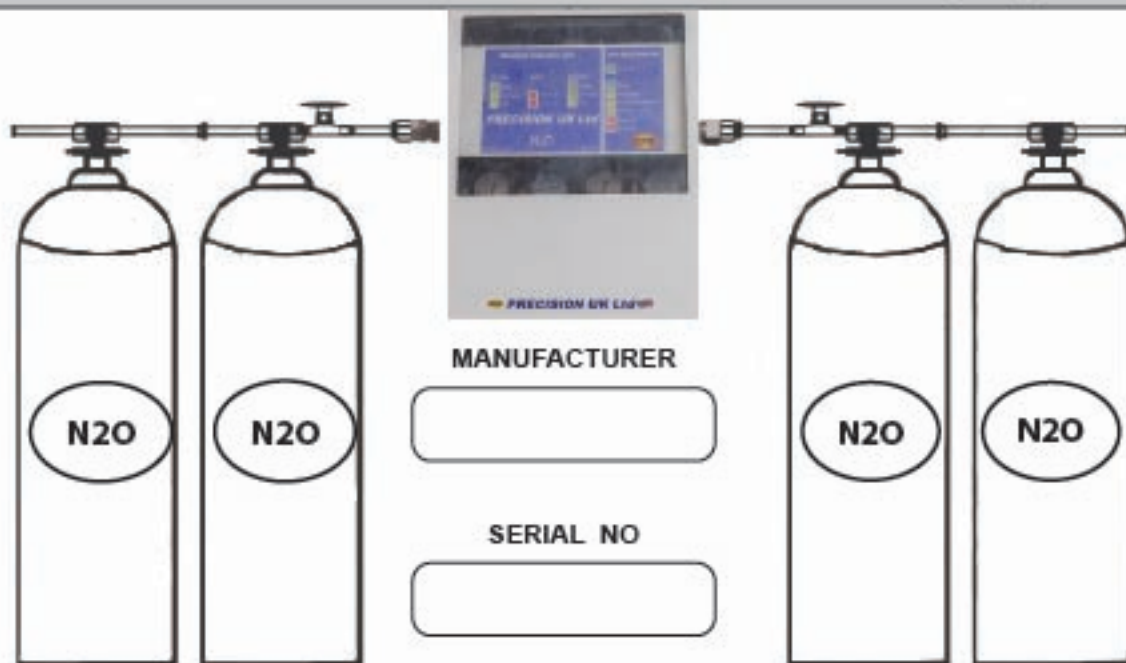
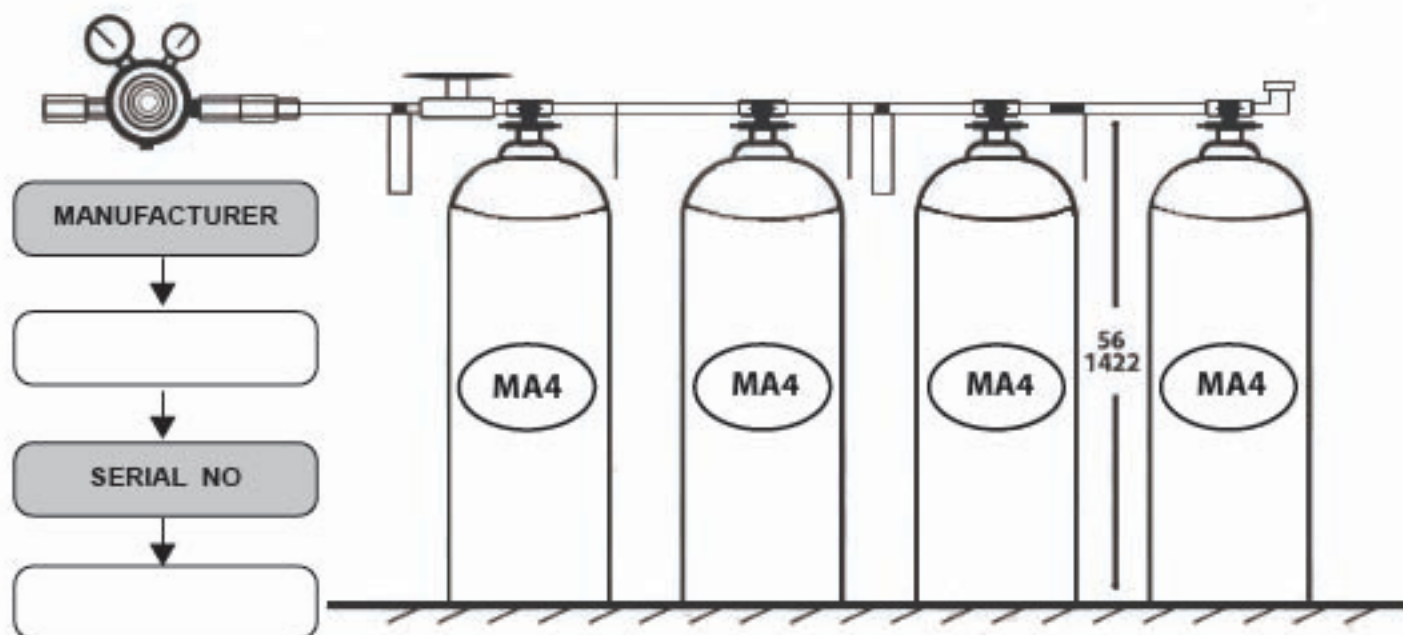
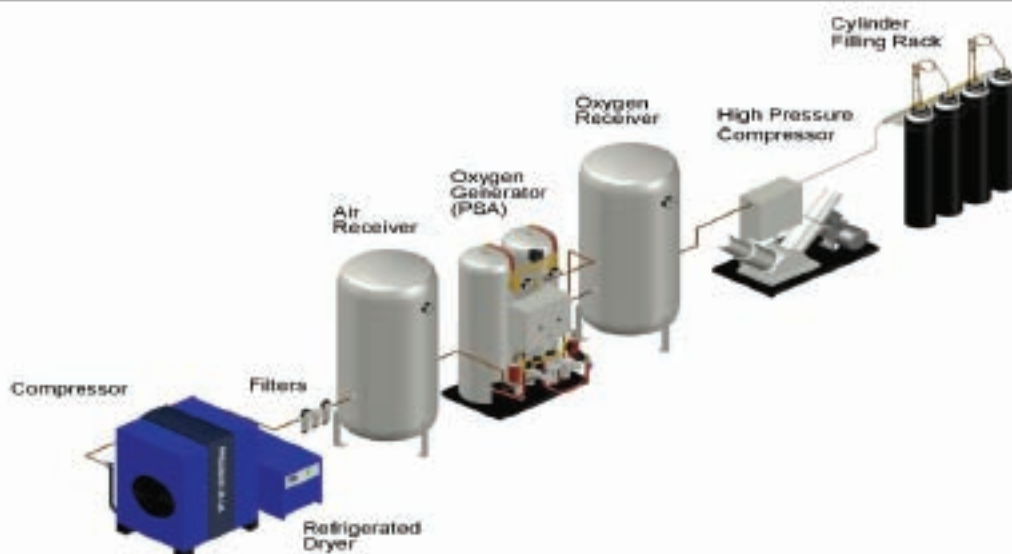
OXYGEN GENERATOR		مولد الأكسجين	
MANUFACTURER			
DATE OF INSTALLATION			
OIL LEVEL			
CONTROL PANEL			
MAX PRESSURE			
MIN PRESSURE			
NOISE			
POWER			

TANK	التانك
MANUFACTURER	
DATE OF INSTALLATION	
LAST PPM DATE	
CAPACITY	
MAX PRESSURE	
DRAIN	
SAFTY VALVE	

FILTERS & SECRATION BOTTLE	الفلاتر
MANUFACTURE	
DATE OF INSTALLATION	
LAST PPM DATE	
LAST DATE OF CHANGE	
CLEAN	
TYPE	
MAX PRESSURE	
MIN PRESSURE	

IGAS ENG NAME:	
CONT ENG NAME:	
M.O.H ENG NAME:	

IGAS ENG SIGNATURE:	
CONT ENG SIGNATURE:	
M.O.H ENG SIGNATURE:	

DATE:

HOSPITAL NAME:

CITY / LOCATION:

CONTRACTOR NAME:

IMPORTANT NOTES

1- You should have stock from oxygen cylinders included regulators , hose , adaptor to connect aneathesia machine , ventilators for Emergency case.

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2- You should have stock from Medical air cylinders included regulators , hose , adaptor to connect aneathesia machine , ventilators for Emergency case.

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3- All dryers should be desicant dryers (chemical dryers) for example if you have refregirant dryers you should chang desicant

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4 - You should have breathing medical air filter.

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5-filter elements should bechanged every 12 months/8000 hours (which ever comes first) activated carbon filter element should be changed every 6months/ 1000hours (whichever

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6-You should change oil every six months or each 3000hours(for vacuum,compressors)

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7- It must be emphasized that the oxygen regulator and the suction regulator should be closed after use to avoid a leakage of oxygen gas and operate the suction pumps unnecessarily.

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8- It is necessary to emphasize the cleanliness of the compressor room and good ventilation in order to enter fresh air or renew the air inside the room, and it is preferable to install air conditioning.

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9- The alarm system should be tested in floors and equipment rooms every 6 months.

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Oxygen
Nitrous Oxygen
Air
Vacuum
Liquid Oxygen Tank
Oxygen Generator

Cairo - Egypt : 21 Hafez Al - Abhar Street
02 / 4872003 02 / 26821616
01271332581 - 01032744799 - 01009084613 - 01015013648
Giza - Egypt : El Nozha Boulak El Dakrour